
# Labofuge<sup>®</sup> 300

# SERVICE MANUAL

P/N 12003230

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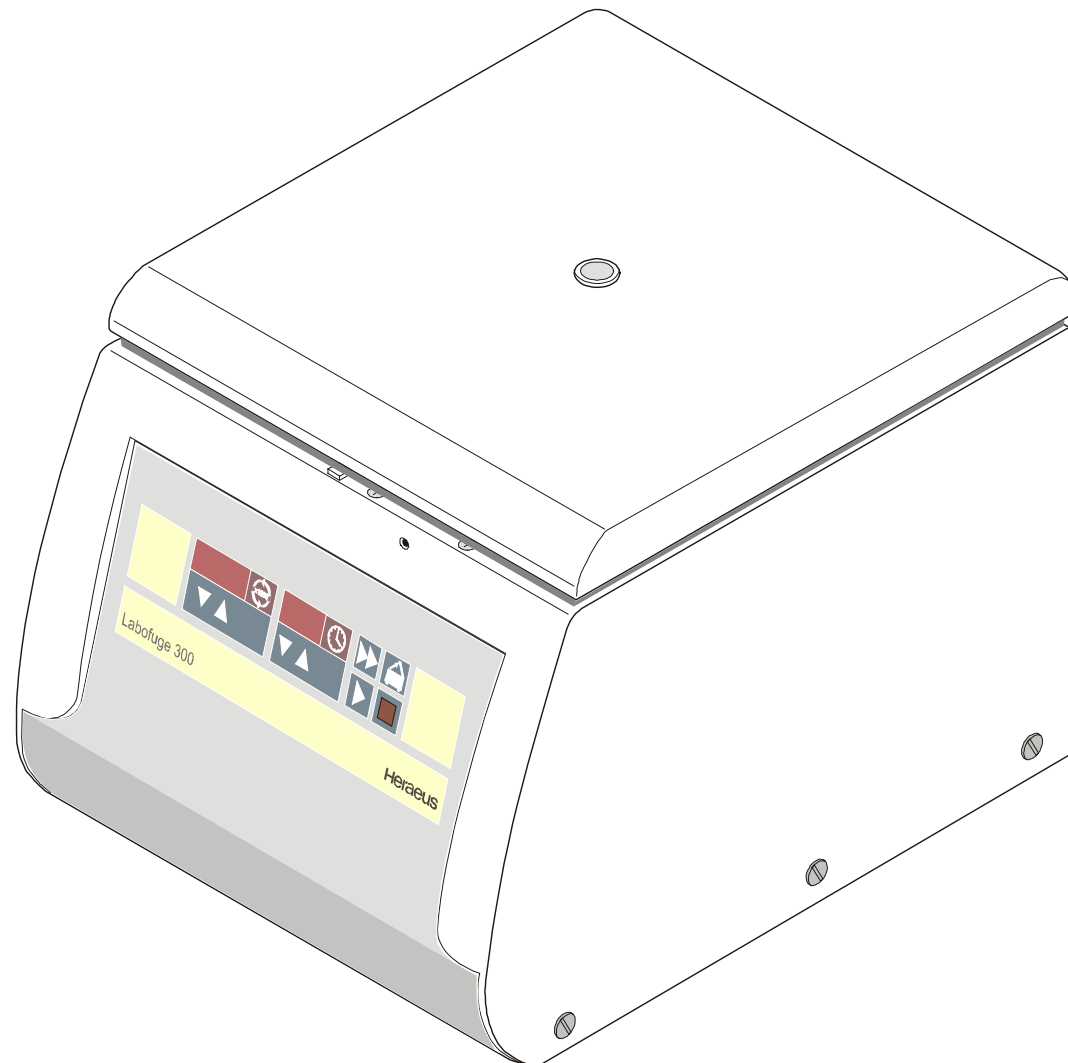**TABLE OF CONTENTS**

Sect.	Title	Doc.- No.	Page
<b>1</b>	<b>OPERATING INSTRUCTIONS</b> (not part of this manual)		
<b>2</b>	<b>SERVICE</b>	Lab300_2	
2.1	Servicing Schedule	" - "	2-1/2
2.2	Trouble Shooting	" - "	2-3/6
2.3	Test Points	" - "	2-7/8
2.4	Imbalance Behavior	" - "	2-9
2.5	Cleaning of Instrument Parts	" - "	2-10
2.6	Electrical Safety Check	" - "	2-11
<b>3</b>	<b>FUNCTIONAL DESCRIPTION</b>	Lab300_3	
3.1	Block Functions	" - "	3-1
3.2	Functions of Main Board	" - "	3-1/2
3.3	Key and Indication Board (Part of Main Board)	" - "	3-3
<b>4</b>	<b>DIAGRAMS</b>	Lab300_4	
4.1	Block Diagram	" - "	4-1
4.2	Wiring Diagrams	" - "	4-2/3
4.3	Wiring Connection Diagrams	" - "	4-4/5
4.4	Main Board 165 (100V) old version - Component Plan	" - "	4-6
4.5	Main Board 165 (100V) new version - Component Plan	" - "	4-7
4.6	Main Board 166 (120V) - Component Plan	" - "	4-8
4.7	Main Board 167 (230V) - Component Plan	" - "	4-9
4.8	Main Board 193 (120V) - Component Plan	" - "	4-10
4.9	Main Board 194 (230V) - Component Plan	" - "	4-11
4.10	Main Board - Wiring Diagram	" - "	4-12/18
<b>5</b>	<b>DISASSEMBLY OF INSTRUMENT PARTS</b>	Lab300_5	
5.1	Housing / Casing Parts	" - "	5-1/2
5.2	Electrical Components	" - "	5-3
5.3	Drive Components	" - "	5-4
5.4	Imbalance Switch	" - "	5-5
<b>6</b>	<b>SPARE PART FIGURES AND LISTS</b>		
6.1	Break Down Drawings	Lab300_6	6-1/6
6.2	Spare Part Lists		
<b>7</b>	<b>Preventive maintenance checklist</b>	Lab300_7	
	<b>Calibration Certificate</b>	" - "	
<b>8</b>	<b>TECHNICAL BULLETINS</b>	Lab300_8	

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Printed in Germany – 3C 05/99 4t Künzel

# LABOFUGE 300 SMALL MEDICAL CENTRIFUGE

## Swing-out value

### Economical and comprehensive

The Labofuge 300 is a small medical centrifuge with swing-out rotor. Included in the price is a rotor with buckets for tubes up to a volume of 15 ml. The Labofuge 300 is ideal for centrifugation of blood collection tubes. It is ideal for medical laboratories and practices with low sample throughput.

### Applications


- centrifugation of blood collection tubes in wards, practices and small hospitals
- preparation of "gel tubes" for virus testing
- separation of cells and crystalline precipitates

### Ideal for virus testing

The Labofuge 300 is particularly suited to diagnostics and therapy of viral diseases (HIV, hepatitis). In the field of virus testing, tubes with a gel barrier ("gel tubes") are increasingly used. Immediately after a blood sample has been taken, these tubes must be spun. With gel tubes, blood plasma is separated from red blood cells by a layer of gel after centrifugation. The gel layer prevents mixing and resultant decomposition of viruses. The tubes can be stored or transported to a laboratory for virus testing.



Centrifuges with fixed angle rotors are unsuitable for gel tubes, according to tube manufacturers.

It goes without saying that ordinary blood sample tubes and glass tubes can also be spun in the Labofuge 300.
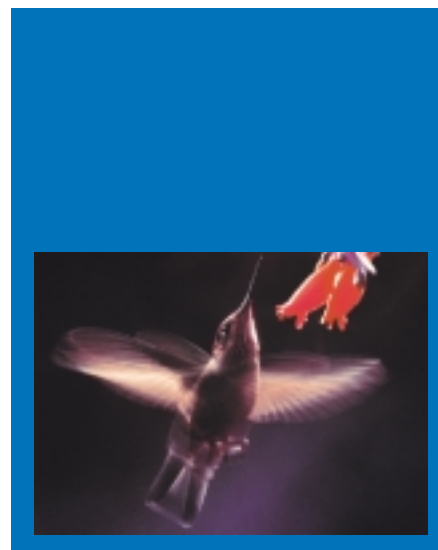
### Benefits

- economical small centrifuge with swing-out rotor
- ideal centrifugation of blood collection tubes
- complete system with 8-place swing-out rotor and buckets
- outstanding quality and reliability as a result of induction drive and microprocessor control
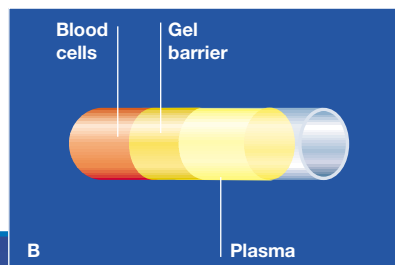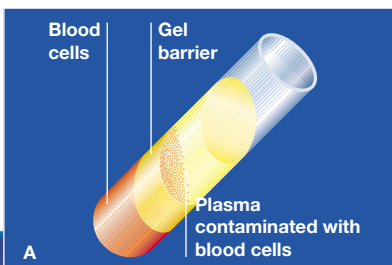
### ◀ Blood collection tubes with gel barrier

- A. Unsuitable for centrifugation in fixed angle rotors.
- B. After centrifugation in swing-out rotor:
  - gel barrier sticks to the tube wall
  - transportation stability is guaranteed
  - plasma has no traces of contamination

- ◀ Swing-out value – with the Labofuge 300 Small Medical Centrifuge

# ORDER NO.

# TECHN. DATA

Type	Equipment	230 V; 50/60 Hz	120 V; 60 Hz
Labofuge 300	incl. 8 x 15 ml swing-out rotor 8 buckets each for 7 ml and 15 ml and distance pads	75003230	75003231

<b>Centrifuge</b>	compact, air cooled table top centrifuge		
<b>Drive</b>	induction drive without carbon brushes		
<b>Control/operation</b>	microprocessor control through Easycontrol I with variable pre-selection and display of speed and time, Quick-run		
<b>Min./max. speed</b>	rpm	300 – 4,000, adjustable in increments of 10	
<b>Max. RCF</b>		2,028	
<b>Noise at max. speed</b>	dB	< 60	
<b>Run time</b>	mins.	0 - 99, continuous operation	
<b>Program memory</b>	the data last entered remains saved		
<b>Construction</b>	galvanised steel chassis with armoured chamber and plastic casing		
<b>Equipment safety</b>	lid lock and interlock		
<b>Dimensions (h/w/d)</b>	mm	315/380/475	
<b>Weight (excluding rotor)</b>	kg	30	
<b>Power input</b>	W	180	
<b>Test standards</b>	manufactured and tested in accordance with EN 61 010-1, EN 61 010-2-020, EN 50 081-1, EN 50 082-1; does not require testing in acc. with UVV VBG 7z		

Rotor/Order no.		Swing-out rotor /76003265 <sup>1)</sup>	35° fixed angle rotor/75003760
<b>Max. speed</b>	rpm	3,600	4,000
<b>Max. RCF</b>		2,028	1,915
<b>Max. capacity</b>	ml	8 x 15	12 x 15
<b>Accel./braking time</b>	sec.	28/29	30/32
<b>Min./max. radius</b>	cm	5.7/14.0	5.1/10.7

<sup>1)</sup> incl. in the price

	Max. tube size		Tubes	Colour	Order no.
	diam.(mm)	length (mm)	per rotor	buckets	
<b>Accessories for swing-out rotor 76003265:</b>					
<b>15 ml bucket</b>	16.5	113	8	red	–
<b>7 ml bucket</b>	13	113	8	yellow	–
<b>Distance pad</b>	–	90	–	–	76003266
<b>Accessories for fixed angle rotor 75003760:</b>					
<b>7 ml adaptor</b>	13	113	12	white	75003761
<b>10 ml pad</b>	16.5	117	6	–	75003762
<b>Protective sleeves for glass</b>	16	100	12	–	75003763

## Your contact

**SORVALL®**  
**Heraeus**

<b>Austria</b>	Kendro Laboratory Products GmbH · Wien · Tel. +43 (1) 801 40-0 · Fax +43 (1) 801 40 40 · office@heraeus.at
<b>France</b>	Sorvall (France) S.A. · Courtaboeuf Cédex · Tel. +33 (1) 69 18 77 77 · Fax +33 (1) 60 92 00 34
<b>Germany</b>	Kendro Laboratory Products GmbH · Hanau · Tel. +49 (61 81) 35-300 · Fax: +49 (61 81) 35-59 44 · info@heraeus-instruments.de
<b>India</b>	Heraeus Instruments India Pvt. Ltd. · New Delhi · Tel. +91 (11) 618 48 40 · Fax +91 (11) 618 53 97
<b>Italy</b>	AHSI S.p.A. · Cavenago Brianza · Tel. +39 (2) 95 08 11 · Fax +39 (2) 95 08 12 77 · ahsidue@tin.it
<b>Poland</b>	Heraeus Sp. z o.o. · Warszawa · Tel. +48 (22) 663 43 23 · Fax +48 (22) 663 43 25 · heraeus.waw@medianet.pl
<b>PR China</b>	Heraeus Instruments Shanghai · Shanghai · Tel. +86 (21) 6417 1565 · Fax +86 (21) 6240 6044 · higchina@guomai.sh.cn
<b>Sweden</b>	Kendro Laboratory Products AB · Upplands Väsby · Tel. +46 (8) 59 07 21 90 · Fax +46 (8) 59 03 16 00 · info@kendro.se
<b>Switzerland</b>	Kendro Laboratory Products AG · Zürich · Tel. +41 (1) 454 12 12 · Fax +41 (1) 454 12 99 · heraeus-ag@swissonline.ch Kendro Laboratory Products SA · Carouge-Genève · Tel. +41 (22) 343 21 67 · Fax +41 (22) 342 38 31 · heraeus-sa@swissonline.ch
<b>USA</b>	Kendro Laboratory Products L.P. · Newtown · Tel. +1 (203) 270-20 80 · Fax. +1 (203) 270-21 66 · info@sorvall.com
<b>U.K./Ireland</b>	Kendro Laboratory Products Limited · Essex · Tel. +44 (1277) 23 15 11 · Fax +44 (1277) 26 18 56 · sales@heraeus.co.uk
<b>All other International Sales in Europe, Middle East, Africa</b>	Kendro Laboratory Products International Sales · Hanau · Germany · Tel. +49 (61 81) 36-300 · Fax +49 (61 81) 35 59 44
<b>Canada, Asia Pacific, Japan, Latin America</b>	Kendro Laboratory Products International Sales · Newtown · USA · Tel. +1 (203) 270-20 80 · Fax. +1 (203) 270-21 66
<b>Internet</b>	<a href="http://www.heraeus-instruments.com">http://www.heraeus-instruments.com</a> · <a href="http://www.sorvall.com">http://www.sorvall.com</a>

**Kendro**  
Laboratory Products

**Kendro Laboratory Products – evolved out of the merger of Heraeus Instruments and Sorvall**

Should your country not be included in the address list, please contact Kendro Laboratory Products International Sales in Germany or USA.

Subject to change without notice. product design wolf raimann wiesbaden

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**2.1 Servicing Schedule (yearly procedure recommended)****2.1.1** Maintenance Routine without Dismantling the Centrifuge**2.1.1.1 Electrical Installation and Safety**

- Switch OFF the centrifuge and disconnect the unit from power, check voltage supply and mains fusing (16 Amps, slow blow characteristic)
- Check condition of plug and wall socket - (let) replace defective parts
- Check cord condition and fixing / connection - replace or refit it
- Check condition of instrument socket and replace it in case of bad contacts

**2.1.1.2 Location and Mechanical Installation**

- Check the base (ground, table, lorry with lockable wheels etc.) For resonance-free and stable conditions
- Check for a well ventilated place and sufficient distances to walls or adjacent equipment, without exposure to direct sunlight
- Check the leveling of instrument (use a spirit level)

**2.1.1.3 Lid Locking Mechanism and Safety Device**

- Connect the centrifuge to power and switch ON
- Check for easy lid closing and self-acting lid opening - if in disorder, readjust lid's swivel hinge and/or lock assembly
- Check the central rubber gasket for lid sealing and replace it, if damaged
- For checking the electrical safety circuit: start the centrifuge, let it shortly run and stop it, the lid must not be unlocked by the microprocessor until the “end” message is shown on the display - if safety circuit is out of function, replace the main board

**2.1.1.4 Cleanliness of Spin Chamber and Motor Casing**

- Open the lid and remove the rotor (for loosening turn socket wrench in arrow)
- Direction - see rim of rotor chamber)
- Clean the spin chamber with a dry and absorbent cloth (remove all dust and moisture - see also 2.5 Cleaning)
- Check the cleanliness of the motor casing and take care of the annular slot around the motor shaft: **penetrating fluids can damage the upper motor bearing**, remove fluids with an injector and/or absorbent paper

**2.1.1.5 Rotor and Accessories Condition and Sealing**

- Check the condition of rotors and accessory parts (especially all supporting or stressed partitions): the rotor and/or accessory parts must not be used any longer, if there are visible traces of mechanical damage or rust
- Check the condition of rotor and/or accessory sealing and replace them in case of malfunction

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## Servicing Schedule

## 2.1.1.6 Rotor Fixing and Motor Shaft

- Check the trouble-free condition of the locking nut and replace it in case of malfunction
- Check the condition of the drive motor shaft: the centrifuge must not be used any longer, if the drive shaft is damaged (bend, thread is worn out, horizontal grooves etc.)

## 2.1.1.7 Temperature Level

- Check the air inlet underneath the lid and under the bottom plate for free ventilation, insufficient air flow will lead to temperature rise of rotor, motor and electronic parts

## 2.1.1.8 Imbalance Behaviour

- Install available and empty rotors and check the imbalance behaviour with rotor dependant cut off and run through weights (see 2.4) and replace worn out motor rubber mounts, motor or faulty main board

## 2.1.2 Maintenance Routine after Dismantling the Centrifuge Casing

## 2.1.2.1 Motor Supporting Elements

- Check the supporting and damping elements of the drive motor and replace them in case of increased rubber abrasion or abundance of imbalance but at least every 3 years

## 2.1.2.2 Braking Circuit

- Check the function of the braking circuit (warming up of brake resistor, even and noiseless brake effect) and replace defective parts in case of malfunction

## 2.1.2.3 Lead and Screwing Connection

- Check the terminal and plug connections of all leads and on all boards and electrical components, tighten all loosen screwing connections, refit or replace defective parts
- Check the screwing connections of all boards, mechanical and electrical components and re-tighten them if necessary (use screw locking lac for motor mounts and lid lock assembly)

## 2.1.2.4 Protection Earth Core and Grounding Connections

- Check the protection earth core for continuity and all grounding plug connectors
- Check isolation resistance and accessible current (see 2.5)

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## 2.2 Trouble Shooting

Error Indication	Error Cause	Possible Error Source	Corrective Procedure	
Displays remain dark	No mains voltage supply	Mains fuse or circuit breaker failed	Check fuse or circuit breaker, replace or switch on again	
		Defective mains cord or switch or instrument socket	Check instrument cord, switch and socket, replace defective parts	
		Defective unit fuse or fuses on main board	Replace it, if fuse blows again, disconnect electrical parts, otherwise search other faults	
	No low voltage supply for indication board	Faulty connection from CPU to indication board	Check connections on CPU, indicat. Board and connecting leads, replace defective parts	
Faulty indication or CPU board		Replace main board completely		
All display elements are shortly illuminated	CPU starts the program again (reset), may be caused by emi	Interruption of running program	NV-RAM out of socket or not correctly placed	Insert the valid NV-RAM and push it correctly into socket
		Reduced voltage supply (<10%)	If the voltage drops often, use a voltage stabilizer	
Centrifuge makes heavy vibrations	Imbalance run	Bad or missing ground connection	Check all ground connections and screwing of all boards	
		Rotor not symmetrically loaded	Open lid, check rotor loading, close lid again and restart	
		Base is not sturdy enough and comes into vibrations	Change or reinforce the base (table, lorry with lockable wheels, etc.)	
		Centrifuge is not correctly levelled	Level the centrifuge correctly, use a spirit level on top of unit	
		Rotor itself has imbalance (mechanical cause)	Rotor must no longer be used, send back to Kendro for inspection	
Drive shaft bended or rotor fixing is damaged	Centrifuge must no longer be used, replace nut and/or motor			

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Trouble Shooting

Error Indication	Error Cause	Possible Error Source	Corrective Procedure
Drive makes noises-no good separation result	Mechanics	Wear out of motor rubber mount	Replace motor rubber mounts (at least every three years)
		Motor bearing	Change motor completely
	Electrical	Defective terminal connecton, faulty lead or motor winding	Check voltage on motor terminal and winding resistances -see test points on boards
		Defective driving	Replace complete main board
Lid cannot be opened by key pressure at standstill	Lid coil is not or not sufficiently supplied with voltage	Missing mains voltage	Remedy see above, manual opening only at standstill
		PTC resistor has released	After a waiting time of 1-2 minutes press key again
		Faulty driving or triac circuit	Replace the complete main board
	Faulty lid coil	Faulty winding of coil	Replace complete lid lock
	Lid is not correctly locked	Lid bolt is jamming	Push lid into lock and press the key again
Lid is deformed or disadjusted		Readjust the lid centrally	
„Lid“ message appears in speed display	Lid was opened manually during run	Forbidden intervention emergency opening device must only be used at standstill	Close lid immediately, turn power off/on, wait for termination of <b>br</b> phase until <b>end</b> message appears
		Defective micro switch or leads or connectors to micro switch are interrupted	Check leads and connectors to micro switch, in case of a faulty micro switch, replace lid lock device completely
	Protection circuit (15V) for lid and motor temperature control was interrupted during run	Motor over-temperature switch has tripped (no sufficient air flow or motor has run on 2 phases only)	Let motor cool down, then check temperature switch and leads with Ohmmeter
			Check air cooling, circulation
			Check solderless connections XM1-3 and motor windings, replace defective parts

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Trouble Shooting

Error Indication	Error Cause	Possible Error Source	Corrective Procedure
„OP“ message in time display	15V supply circuit is interrupted at standstill	Micro switch	Check leads and micro switch
		Motor over-temperature switch	Check motor over-temp. Switch and leads, replace defective parts
„br“ message appears in speed display	Rotor comes to standstill without braking force	Short interruption of mains supply	Wait for rotor standstill (appr. 75 seconds) and re-start
		Bad contact of instrument supply line	Check plug contacts of supply line and replace defective parts if necessary
„bAL“ message appears in speed display	Imbalance run	Rotor not symmetrically loaded	Open lid, check rotor loading, close lid again and restart
		Base is not sturdy enough and comes into vibrations	Change or reinforce the base (table, lorry with lockable wheels, etc.)
		Centrifuge drive is not correctly leveled	Level the unit correctly
		Rotor itself has imbalance	Rotor must no longer be used, send back to kendro
		Drive shaft	Replace the motor
		Signal fault	Imbalance switch
	Interruption at the wiring		Check wiring to the imbalance switch and replace if necessary
	Circuits on the main board		Replace the main board
	Imbalance but no „bAL“ message	No imbalance cut-off	Imbalance switch
Circuits on the main board			Replace the main board
E- „0“ message appears in speed and time field	Rotor didn't turn	Rotor is jammed	Check for easy rotor movement, remove any jamming objects
		Motor is jammed	Replace motor

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**Trouble Shooting**

Error Indication	Error Cause	Possible Error Source	Corrective Procedure
“E-“ “0“ message appears in speed and time field	Motor didn't start	Connections inter drive and main board	Check terminal and lead connections, replace faulty parts
		Defective drive	Check resistance of motor windings, replace faulty parts
		Faulty main board	Remove main board completely and replace it
	Motor didn't accelerate	False rotation direction	Change 2 motor connection lines XM
“E-“ “8“ message	Over-voltage of intermediate circuit	Defective leads or brake resistor	Replace main board completely
		Driving and/or braking circuits	Replace main board completely
“E-“ “10“ message	Checksum error of NV-RAM	NV-RAM is not initialized or false	Check NV-RAM and socket, insert the correct NV-RAM
“E-“ “11“ message	Disturbed data transfer from NV-RAM	Sporadic: strong electromagnetic disturbance	If error returned very often, screen the disturbance source or operate the centrifuge on another place
		Continuously: error on main board	Remove main board completely and replace it

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### 2.3 Test Points

Test Points	Unit value	Conditions
Mains termin. <b>XN</b> unit's resistance	120/230V 1360Ω	all given values refer either to 100V (Japan), 120V (US) or 230V 50/60Hz, ( $\pm 10\%$ )
Terminal <b>XM</b> Motor voltage - Values must be the same between all 3 motor leads	115V AC 175V AC 215V AC	in each case measured inter 2 motor leads after reaching the selected speed 1200min <sup>-1</sup> , swing-out rotor #76003265 2400min <sup>-1</sup> 3600min <sup>-1</sup>
	95V AC 180V AC 225V AC	1000min <sup>-1</sup> , fix angle rotor PP #75003760 2500min <sup>-1</sup> 4000min <sup>-1</sup>
Motor current <b>I<sub>M</sub></b> Cable 1, 2, 3	3x1,5A	soft iron or digital effective measuring instr. maximum during acceleration (#3265, 30s)
	3x1,2A 3x1,1A	rotor #3265, set speed = 3600min <sup>-1</sup> , 180W rotor #3760, set speed = 4000min <sup>-1</sup> , 130W
	3x1,0A	maximum during braking phase (#3265, 30s)
Intermediate voltage <b>U<sub>D</sub></b> (brake resistor -> current control resistor R44)	320V DC 320V DC	at standstill max. during braking phase
	308V DC 308V DC	rotor #3265, set speed = 3600min <sup>-1</sup> rotor #3760, set speed = 4000min <sup>-1</sup>
	305V DC	shortly at maximum acceleration
Current <b>I<sub>D</sub></b> of intermediate circuit, (measured as voltage drop across control resistor R44) <b>U<sub>M</sub> = I<sub>D</sub> x R<sub>M</sub></b>	25mV DC 60mV DC 100mV DC	1200min <sup>-1</sup> , swing-out rotor #3265 2400min <sup>-1</sup> 3600min <sup>-1</sup>
	20mV DC 50mV DC 70mV DC	1000min <sup>-1</sup> , fix angle rotor PP #3760 2500min <sup>-1</sup> 4000min <sup>-1</sup>
	125mV DC	shortly at maximum acceleration (#3265)
Motor winding resistance 20°C -insulation value	3x20Ω > 10MΩ	switch OFF unit, pull off motor plugs, measure inter 1 -> 2, 2 -> 3, 3 -> 1 resistance inter each phase and motor casing
Brake resistor R51 soldered	220Ω	switch OFF unit, resistance at 20°C
Lid solenoid Connectors <b>XC1 XC2</b>	210Ω 40Ω	Old version without temperature switch resistance at 20°C 230V version, p/n 20901100 120V version, p/n 20901101

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**Test Points**

Test Points	Unit value	Conditions
Lid solenoid Plug <b>XC</b>	80Ω 23Ω	New version with temperature switch resistance at 20°C 230V version, p/n 20901130 120V version, p/n 20901131
Lid micro switch Plug <b>XA</b>	15V DC	voltage drop by open lid
Motor over-temp. Connectors <b>XB</b>	15V DC	voltage drop by closed lid and one disconnected lead
Imbalance switch Plug <b>XE1 XE2</b>	5V DC	voltage drop if imbalance switch in neutral position

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## 2.4 Imbalance Behaviour

### 2.4.1 Meaning of Imbalance Switch

**The imbalance switch has the task to switch off the centrifuge by incorrect loading of the rotor before reaching the normal speed.**

#### NOTE!

- Before starting with the re-adjustment procedure check the greasing of rotor trunnions and cleanliness of bucket grooves. Opposing buckets and racks must be identical and equal in weight
- Check the motor's rubber mounts and the leveling of the drive with installed rotor (use a spirit level)

### 2.4.2 Adjusting Procedure

For adjusting the imbalance switch a 1,5mm gauge is required. For the gauge a washer is used, which is located near the empty space at the ground plate under the gas lid stay.

#### Sequence of operations:

- Loosen the two screws at the ground of the mounting bracket, so that the mounting bracket can be move
- Put the gauge between the mass ring and the roller at the micro switch, take care that the gauge is at right position
- Move the bracket away from the mass ring, then move it backwards to the mass ring, so that the switch point of the micro switch will be reached. If the motor will be moved, the micro switch will switch lasting to ON or OFF position
- Tighten the screws at the mounting bracket and secure them with screw locking lack

### 2.4.3 Test Run Performance

- **Install** in table available **rotor(s)** in unloaded condition

Rotor	Run through weight	Cut off weight
75003265	approx. 8g	approx. 15g

- Insert rotor dependent **cut-off weights** and position them at angles of 90° to each other, the centrifuge must stop 4 times indicating “bAL” message - in case of false behaviour see trouble shooting: no “bAL” indication
- Insert the **admissible imbalance weights** in the same manner, the centrifuge must run through 4 times to maximum rotor speed

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## 2.5 Cleaning of Instrument Parts

### **ATTENTION - WARNING!**

**The electrical and electronic components must not be cleaned with moist detergents!**

For Cleaning the centrifuge housing or its accessories see Operating Instructions section 1 (maintenance and care).

- **Electronic components**

Clean dusty components carefully with a dry and soft brush and remove loose dust with a vacuum cleaner

- **Vent holes**

Remove dirt from vent holes of lid or bottom plate using a brush and vacuum cleaner

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## 2.5 Electrical Safety Check

### ATTENTION!

A final electrical safety check must be performed after each maintenance and/or repair!

- **Resistance check of protective conductor**


The measuring value of the resistance between the mains plug's grounding pin and the grounding conductors of the motor, electronic chassis and the casing must not exceed 200 mΩ.

- **Insulation resistance Check**

Check also the insulation resistance between the poles of the mains plug and the grounding conductor; the resistance value must be more than 2 MΩ.

- **Accessible current measured to EN 61 010**

The accessible current must not exceed 3.5 mA in single fault condition (interrupted protection earth wire)! In accordance with the EN61010, IEC1010 and UL3101 such a fault condition can be reproduced by the following measuring circuit.

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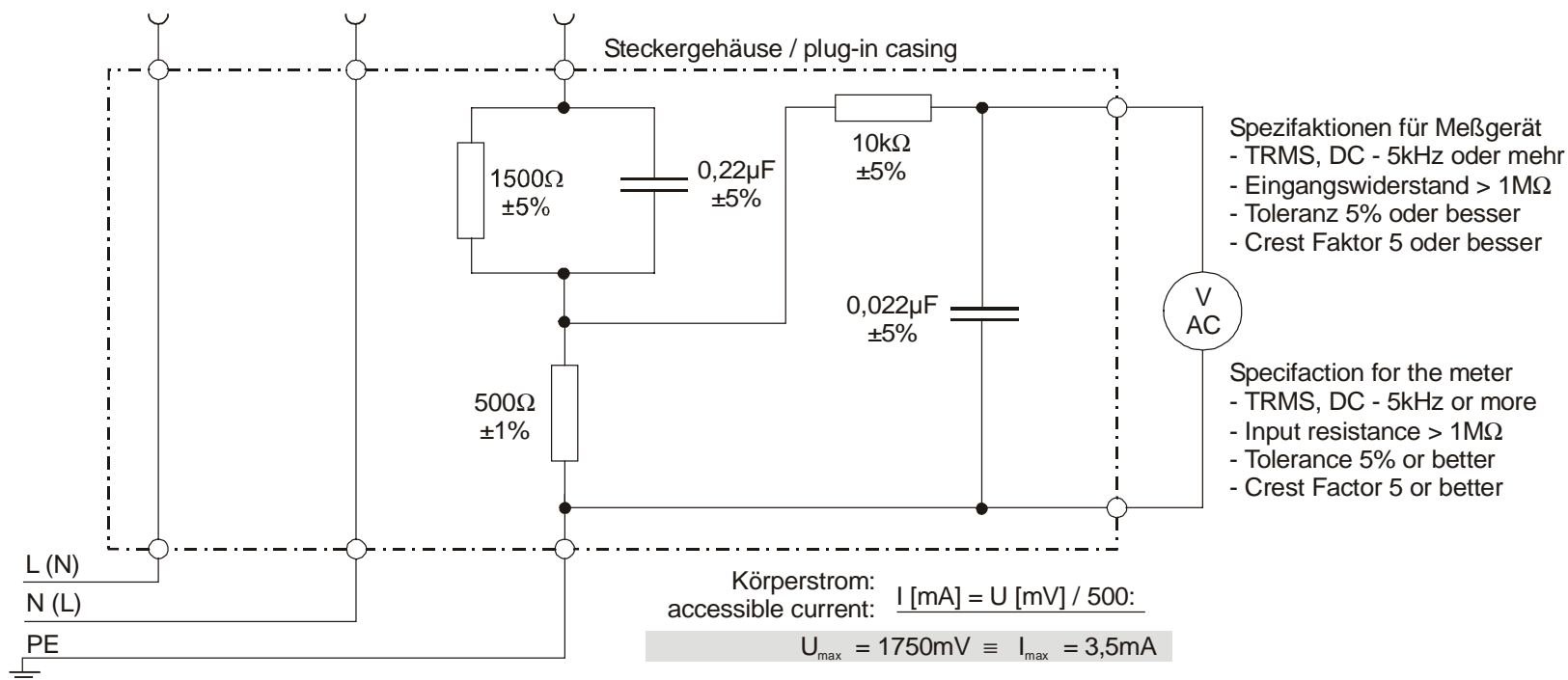### 3.1 Block Functions

The Die Labofuge® 300 is a microprocessor controlled laboratory tabletop centrifuge with induction drive motor and integrated air cooling system.

The unit incorporates following boards and components (see block diagram 4-1):

- Two poles mains switch with integrated instrument plug for mains cable mounted on the unit's back side
- Main board with microprocessor part and power electronics, 2 fuses are serving for board protection
- Key and indication board (EASYCONTROL programming without temp.display)
- 3 phase induction motor with integrated thermal overtemp. Switch (C. O. 140°C)
- Lid lock assembly with solenoid and integrated micro switch (mechanical bolt keeper, magnetic unlocking), mounted under the casing's right side
- Brake resistor (part of main board) is mounted onto the ground plate
- Imbalance switch is mounted onto the ground plate

### 3.2 Main Board Functions

The main board is mounted on the ground plate in front of the drive.

The components on main board are arranged in following groups (see wiring diagram page 4-2 and 4-3).

- Fusing (2 x 4/2 Amps. Slow blow), noise filter in accordance with EN 55011
- Power pack for low voltage supply of microprocessor part (physically separated by transformer) and power electronics (mains potential)
- Triac control circuit for lid solenoid ( DC supply via bridge rectifier)
- Bridge rectifier for DC intermediate circuit supplying brake control path and frequency converter
- Microprocessor part with controller (CPU) and ASIC80 component
- Exchangeable NV-RAM containing specific data
- Power part inclusive physically separated driving stages

#### 3.2.1 Power Pack

The power pack (transformer, bridge rectifier and voltage regulator) generates:

- U1 = 5 V supplies the CPU part on main and indication board  
Reference potential **A1**: connected to protective conductor (GND)
- U2 = 15 V (generated by one way rectifying) supplies lower HGTPs of power Electronics via the protection circuit of closed lid switch and motor overtemperature switch (normal closed)  
Reference potential: **A2! Mains potential!**  
**EXERCISE CAUTION when measuring etc.**

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### Main Board Functions

#### 3.2.2 Intermediate Circuit with Brake Path and Frequency Converter

The DC intermediate circuit serves as an energy store between the AC power input and the transmitted motor performance. The intermediate circuit consists of a bridge rectifier (230V) and reservoir capacitors which are softly charged via NTC resistor, continuously supplied during operation and slowly discharged after power off. The 100/120V variant has a voltage doubling stage (two diodes for charging two serial connected reservoir capacitors).

##### Brake Path

Electrical power can be fed back into the intermediate circuit during motor deceleration (motor acts as generator). This braking power is transformed into heat by the line soldered resistor so that the intermediate circuit voltage does not rise to an excessive level. The brake resistor is switched into the intermediate circuit by a fast switching transistor (pulse-width modulation, synchronized by mains frequency). This FET is voltage dependent controlled by a self-acting stage (closed loop).

##### Frequency Converter

The motor is 3-phased provided with chopped direct voltage blocks. These blocks are variable in frequency and pulse-width modulation, dephased to 120°. These 3 phases are controlled during acceleration, running at set speed and deceleration (e.g. for small speed low frequency and small pulse-width length will be affected) The frequency converter is protected against over-current and over-voltage.

#### 3.2.3 Microcontroller (Central Processing Unit) Part

The software identification numbers of the CPU 80C52 and data storage (NV-RAM) 9346 are displayed in speed and time fields in sequence:

- CPU: 478 xx
- NV-RAM: 1067 xx

The actually indicated version numbers (xx) are displayed in the time field and will be upgraded if program or data changes are necessary. The control program (ROM) is integrated in the controller component. The most important operating parameters (e.g. the maximum speed and the last operator settings) are stored in the non volatile (NV)-RAM. The CPU cooperates with the ASIC 80 component via the 8 bit data bus and some control lines. The ASIC 80 is a specially designed integrated circuit for the driving of Kendro centrifuges with induction motor, containing some protecting functions of the power electronics. The CPU signal to unlock the lid is combined with the reset controller chip, this measure makes sure that the lid cannot be opened in case of a power failure.

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### 3.3 Key and Indication Board (Part of Main Board)

The key and indication board is mounted behind the operating panel. The connection to the main board is done by a 14 polar flat cable. The 7-segment displays, the control LED's and the operating keys are managed by the CPU by multiplex processing. For the programming serves EASYCONTROL (see Operating Instruction Manual).

#### Function of keys and indicating elements

Keys	Name	Function(s)
S1	Lid symbol	To open the lid (only possible at standstill) – else error behavior (PTC is heated up)
S2	Double arrow	Quick run, to start - keep the key pressed
S3	Square symbol	Stop, to terminate the centrifugal run manually
S4	Arrow to right	Start, to set a centrifugal run in rotation
S5	Speed arrow down	To decrease the value of the set speed
S6	Speed arrow up	To increase the value of the set speed
S7	Time arrow down	To decrease the value of the set time
S8	Time arrow up	To increase the value of the set time

Indicator	Name	Function(s)
H1,H2	Time range	Indication of centrifugation time and hold (hd), <ul style="list-style-type: none"> <li>• Right part of software identification numbers of CPU and NV-RAM</li> <li>• Actual version numbers (“xx“)</li> <li>• Lid OPEN condition (“OP“)</li> <li>• Code numbers of error messages</li> </ul>
H3 - H5	Speed range	Indication of speed [x1000 (rpm)] <ul style="list-style-type: none"> <li>• Left part of software identification number</li> <li>• Run finished (“End“)</li> <li>• Error messages like (“lid“, “E-“)</li> </ul>
H6 - H9	Rotation LED's	Indication of rotation after starting till “End“

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The index numbers stated in (xxx) reappear within the break down drawings and spare part list.

## 5.1 Dismantling the Housing

### 5.1.1 Outer Casing (100)

- Open the lid, switch off the unit and **pull out the mains plug**
- Unscrew and remove the rotor
- Unscrew and remove the plastic motor cover (102)
- Disassemble the lid see 5.1.2
- Unscrew 6 flat-head screws (104) for bottom plate attaching and remove them together with plastic washers (103)
- Lift the casing in front and disconnect the wires of the lid lock assembly (be carefully when pulling off the connectors XC1,2 with nose pliers), or unplug the plug XC on units with new lid lock
- Remove the casing completely with lid lock assembly in vertical position across armored chamber and lay on side
- Refit the casing analogously in reverse order

### 5.1.2 Lid

Note! The lid consists of the outer plastic casing (110) and of the screwed-in steel plate (112) with it's bended tabs used for keeping the lid closed and for attaching the hinges (109) and to fix the fork-type joint (114) for gas lid stay.

- Open the lid, switch off the unit and pull out the mains plug
- Knock carefully the joint bolt (115) for upper attaching of gas lid stay (116) out of fork-type joint (114) with its three notches ahead by use of hammer and pin punches (but press the bolt in with non notches ahead by use of pliers)
- Push an adequate screw driver into the bolt's place to keep the gas lid stay in position and to avoid lid's slamming
- Remove the 4 screw caps (107) on unit's back side to get access to the 4 screws (108) for both hinges' attaching and use a 3mm Allen key for unscrewing them (if the casing borings should not be in alignment to the screw heads, then loosen the 6 flat head screws (104) and press the casing in the correct position to remove the screws)
- Take the screw driver out of fork-type joint to separate lid from gas stay and remove the lid from casing
- For adjustment / exchange of lid hinges (109) or the fork-type joint (114) remove the 4 screws (113) for the lid's lower metal plate attaching

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### 5.1.3 Lid lock (140)

Note! Unlocking solenoid and micro switch are components of one single spare part and therefore cannot be ordered separately!

- Dismantle the casing - see 5.1.1
- Take off the lines for solenoid (XC1, 2 / XC) and microswitch (XA) at the main board
- Loosen both counter-sunk screws to take out to both sides the lock's clamping nuts on guide plates
- Remove the screws and take the lid lock device out to the right side guiding it between front plate and armoured chamber
- Open the emergency cord and pull out of it's guide holes
- When reassembling the lock adjust it centrally to bolt keeper's position

### 5.1.4 Armored Chamber

- Dismantle the casing - see 5.1.1
- Remove the 6 hexagon screws for bottom plate attaching by using a small 10mm socket wrench with at least 26cm of length
- Expose connecting lines to indication board from insulation paper (143)
- Lift armoured chamber out of guide slots in bottom plate and lay down on it's right side
- When re-assemble pay attention to the correct seat of the insulation paper!

### 5.1.5 Gas Lid Stay (116)

open the lid, switch off the unit, pull out the mains plug and remove the rotor

- Knock carefully the joint bolt (115) for upper attaching of gas lid stay (116) out of joint block (114) with its three notches ahead by use of hammer and pin punches (but press the bolt in with non notches ahead by use of pliers)
- Push an adequate screw driver into the bolt's place to keep the gas lid stay in position and to avoid lid's slamming
- Lay unit carefully on it's right side and remove both hexagon screws for lower joint block attaching
- Take screw driver out of upper block and pass the gas lid stay together with lower block downwards through opening of bottom plate
- Remove the lower joint bolt as described above
- Install the new gas lid stay analogously in reverse order

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## 5.2 Replacement of Electrical Components

### 5.2.1 Fuses on Main Board

- Dismantle the casing - see 5.1.1
- Check both micro fuses and if necessary replace them by equivalent and undamaged ones (see schematic diagram 4 - 2)
- Perform a test run and search for blowing cause
- Reassemble the device analogously in reverse order

### 5.2.2 NV-RAM (122) on Main Board

- Notice old displayed NV-RAM's version number when powering the unit on (2. No. in sequence)
- Remove the casing (5.1.1), if necessary dismantle armored chamber (5.1.4)
- **ATTENTION - CMOS components! Discharge your body before handling!** Notice correct position of NV-RAM (8 pins) and pull carefully out of socket
- Reinsert the new NV-RAM correctly
- Reassemble the device analogously in reverse order, switch unit back to power and check the displayed NV-RAM's version number

### 5.2.3 Main Board (120) Replacement

ATTENTION! The key and indication board as well as the brake resistor are parts of the main board and cannot be ordered separately!

- Dismantle the armored chamber as described in 5.1.4
- Disconnect all single plug connectors as: mains supply **XL, XN**, lid solenoid **XC1,2**, motor lines **XM1-3** and over-temperature cut-out **XB1,2** by use of nose pliers and remove the plug for lid switch **XA**
- Remove all screws of the main and indication board as well as brake resistor, remove the washers and take out the parts
- Pay attention to the notes in 5.2.2 and remove the inserted but non programmed NV-RAM of the new board
- Take the new NV-RAM (122) out of transport box or re-use the old but **faultless** component and insert it into the socket of the new main board
- Separate indication board from new main board on its' breaking line and mount these parts with all washers, screws and/or nuts on its' correct positions as well as the brake resistor
- **Do not mix up disconnected cables during re-connection!**
- Reassemble the device in reverse order and perform a test run, **making sure the drive turns in the right direction** (see imprinted arrow on rim of rotor chamber)!

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### 5.3 Replacement of Drive Components

#### 5.3.1 Disassembly of Drive Motor (130)

- Remove the casing (see 5.1.1)
- Disconnect the leads for motor **XM1-3**, over-temperature switch **XB1,2** from main board by use of flat pliers and unscrew the grounding wire **PE** from motor casing
- Unscrew the 3 hexagon screws for motor attaching on mass ring and remove them together with the lock washers
- Lift the motor out of mass ring and take it out
- Reassemble the motor analogously in reverse order, replace the lock washers, fasten the 3 nuts and fix them with screw locking lack (138). Perform a test run. Check the sense of rotation (see arrow direction on rim of rotor chamber) and mix up 2 leads on terminal **XM1-3** in case of disorder

#### 5.3.2 Motor supports (132)

**All 3 rubber supports (132) have to be replaced at the same time and at least every three years! By the way the all double-sided sandpaper discs (133) must be exchanged, too.**

- Open the lid, switch off the unit and **pull out the mains plug**
- Unscrew and remove the rotor
- Unscrew and remove the plastic motor cover (102)
- Unscrew the mass ring from rubber supports, remove attaching nuts and lock washers
- Hold motor and mass ring on top of rubber supports with one hand and lay the unit down onto the right side (e.g. on a towel)
- Pull motor together with mass ring lateral out of rubber supports and lay down carefully inside rotor chamber onto a cloth
- Unscrew the rubber supports underneath bottom plate and remove them
- Install the new rubber mounts and sandpaper disks, push the motor back on top of rubber mounts and hold in position when putting the unit back on it's feet
- Replace lock washers, tighten the nuts evenly (10 Nm, when using a torque key) and secure them with screw locking lack (138)
- Reassemble the device in reverse order. At last perform a test run

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## 5.4 Imbalance Switch

### 5.4.1 Replace the Imbalance Switch

- Remove the casing (see 5.1.1)
- Dismantle the armored chamber (see 5.1.4)
- Remove plug **XE** at the main board
- The imbalance switch consists out of a special micro switch with wiring (170) and the mounting bracket. For replacement unscrew the 2 screws at the switch
- Replace the switch (170), tighten the screws and secure them with screw locking lack (138), replace the plug **XE**
- After each modification at the drive system (replacement of parts) and imbalance switch several test runs have to be done to make sure that the unit allows maximum possible imbalance for the customer – but must cut off not allowed imbalance

### 5.4.2 Imbalance Switch adjustment


- For adjusting the imbalance switch a 1,5mm gauge is required. For the gauge a washer is used, which is located near the empty space at the ground plate under the gas lid stay
- Sequence of operations: Loosen the two screws at the ground of the mounting bracket, so that the mounting bracket can be move
- Put the gauge between the mass ring and the roller at the micro switch, take care that the gauge is at right position
- Move the bracket away from the mass ring, then move it backwards to the mass ring, so that the switch point of the micro switch will be reached. If the motor will be moved, the micro switch will switch lasting to ON or OFF position
- Tighten the screws at the mounting bracket and secure them with screw locking lack

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
- Insert rotor dependent **cut-off weights** and position them at angles of 90° to each other, the centrifuge must stop 4 times indicating “bAL“ message - in case of false behaviour see trouble shooting: no “bAL“ indication
- Insert the **admissible imbalance weights** in the same manner, the centrifuge must run through 4 times to maximum rotor speed

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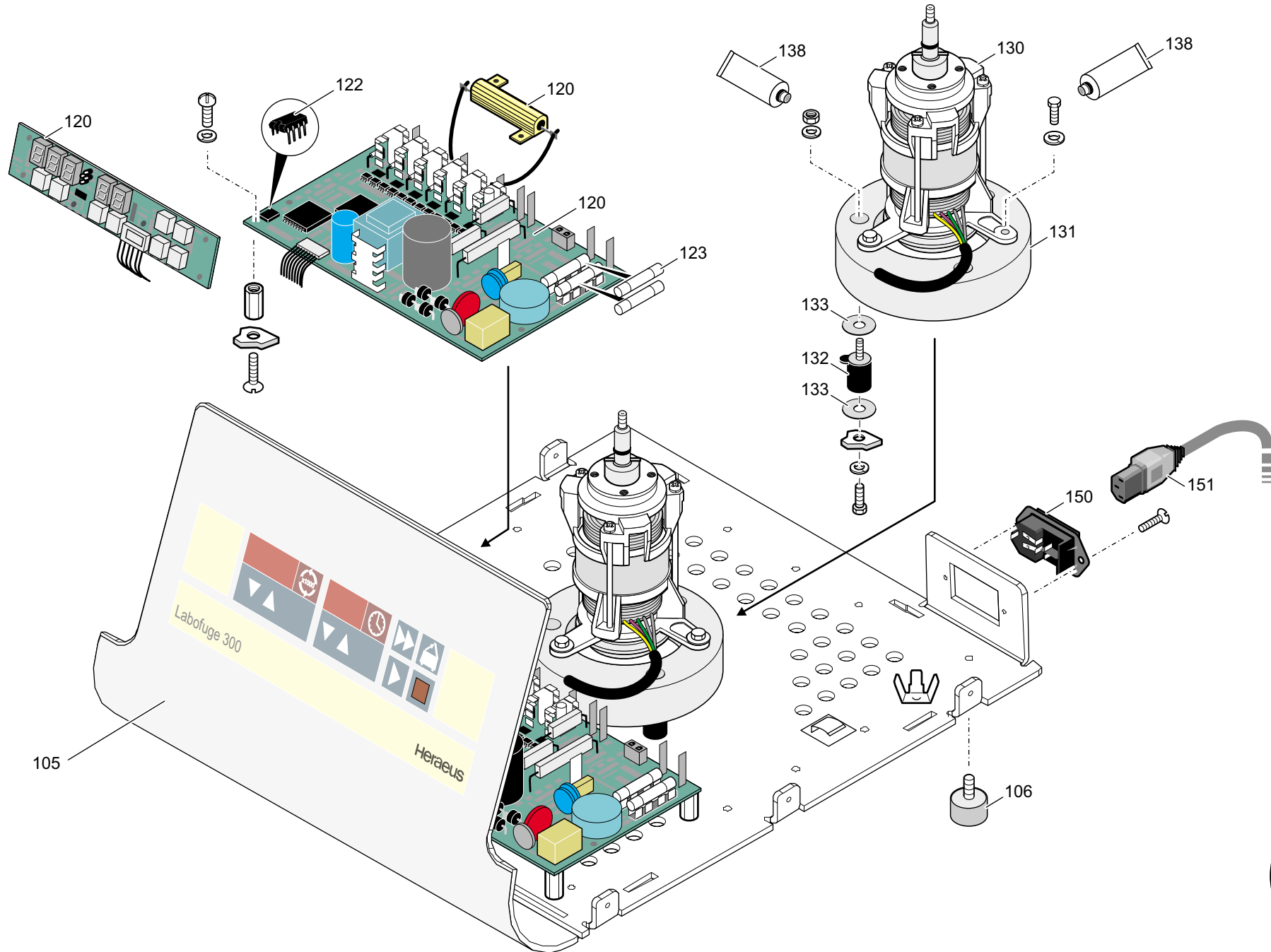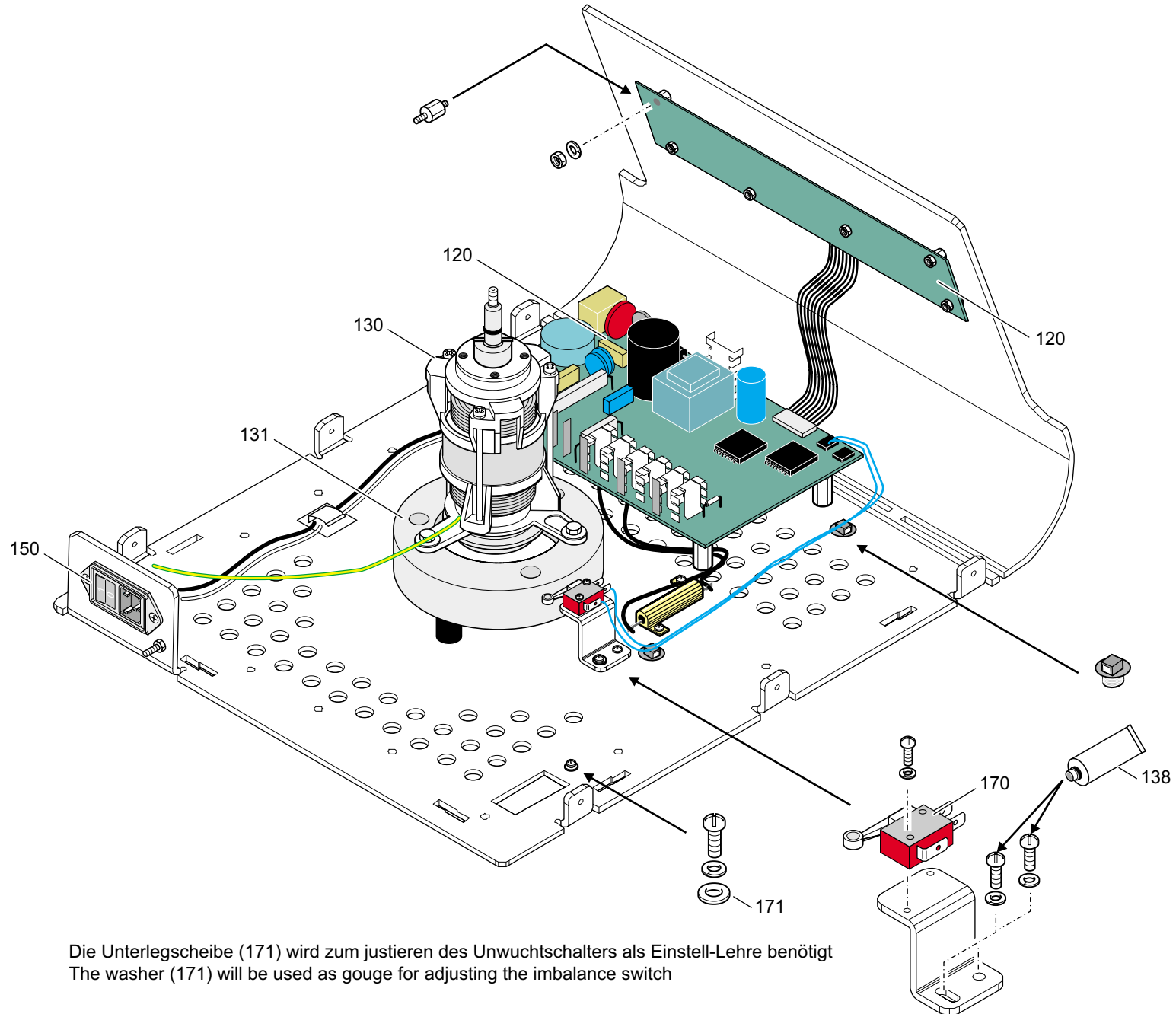# Calibration Certification Procedure

## Labofuge 300

Revised – 13/05/02

### 1. General

Successfully completing the Preventive Maintenance Checklist and Calibration Certification Procedure as outlined will ensure that the instrument is properly maintained and calibrated to Kendro specifications.

### 2. Equipment Required

- Digital Multimeter
- Stopwatch
- Photo Tachometer
- Safety Tester
- Labofuge 300 Service Manual
- Centrifuge Operating Manual

**Note:** Check calibration due date on all test equipment before starting certification.

### 3. Documentation Required

- Preventive Maintenance Procedure and Checklist (Service Manual)
- Calibration Certification Procedure and Form
- Certification Label
- Certificate of calibration for each piece of test equipment used to perform the certification
- Training Certificate

#### Preventive Maintenance Check

- Perform Preventive Maintenance checks as outlined in the Kendro Service Manual to ensure the instrument is in good working order without performing any calibrations.
- Complete the Preventive Maintenance Checklist. (calibration data will be filled in after Certification)

#### Calibration Certification

- **SPEED**

Install a rotor, set the speed control to 3000 rpm and start the centrifuge. When the speed has stabilized, measure the actual speed with an external tach (photo tach on multimeter where speed = frequency[no conversion]). Record measured speed and displayed speed in the speed section of the certification form. Set the speed control to max. rpm and after stabilization, record measured speed and displayed speed in the speed section of the certification form. Determine if the speeds recorded meet the specifications on the form. If the specifications are met, proceed to the next section. A recalibration of the instrument is not possible.

# **Calibration Certification Procedure**

## **Labofuge 300**

- **TIME**

Using a stopwatch, check the timer by setting a 10 minute run. Start the centrifuge and measure from the time the start button is pushed until the timer switches to stop. Record the measured time in the time section of the form. Determine if the time value recorded meets the specifications in the form. If the specifications are met, proceed to the next section. A recalibration of the instrument is not possible. Initial the time section.

- **SAFETY TEST**

Using safety tester and connect the centrifuge to it. Measure Insulation resistance.. Select to Earth Conductor resistance and touch with probe all mechanical part of the centrifuge. Resistance must be <0.2 ohms at all points otherwise check for good grounding. Measure accessible current to procedure in the service manual.

### **Documentation**

After filling in the certification information complete the following documentation:

- Log test equipment data on calibration certification form.
- Review P.M.checklist and calibration certification form to ensure that all data and required information has been recorded properly.(including your signature and date)
- Fill out a certification label and place it on the right side of the centrifuge where visible or at a location selected by the customer.
- Ask the customer to sign the calibration certification form and place a copy in the customer's Certification Record Book for that centrifuge.
- Place copies of the P.M.checklist and Field Service Report in the same Record Book.
- Provide copies of your training certificate and test equipment calibration certificates if requested by the customer

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Kendro service Information

# Heraeus Centrifuges

## PREVENTIVE MAINTENANCE CHECKLIST

Labofuge 300

Type  
S/N

### Rotor Inspection

- o Check for corrosion and wear
- o Inspect rotor lid locking screw
- Serial numbers:  
75003265.....  
75003760.....  
.....  
.....  
.....  
.....  
.....  
.....

### Performance Checks

- o Acceleration time \_\_\_\_\_ seconds
- o instrument top speed \_\_\_\_\_ rpm
- o Deceleration time \_\_\_\_\_ seconds
- o Time (10 minute set) \_\_\_\_\_ actual min:sec

### Imbalance

- o Centrifuge stable and level
- o Centrifuge shuts off w/proper weight

### Preliminary Checks

- o Rotor chamber (defrost)
- o Inst. environment  
(level & clearance)

### Safety test

- o Accessible leakage current <3500µA
- o Insulation resistance >10 Megaohm
- o Earth Conductor Resistance <0.2 Ohm

### Pre-run Checks

- o Ground continuity
- o Line voltage \_\_\_\_\_ V AC
- o Display segments on power-up
- o Door interlock and latch
- o Drive stop & Door unlatch (50 rpm)

### Mechanical

- o Motor bearings
- o Microprocessor Initialization (if required)

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Standard Test			Rotor Used:			
Set Speed (rpm)	Tolerance (whichever is greater)	Displayed Speed (rpm)	Measured Speed (rpm)	Meets Spec. (y/n)	Measured Speed after Cal.(rpm)	Passed (Initials)
2000	120 rpm					
4000	120 rpm					
Custom Test			Rotor Used:			
	120 rpm					
	120 rpm					

#### TIME / ZEIT

Standard Test		Rotor Used:			
Timer Set (min.)	Tolerance (sec.)	Measured Time (m/s)	Meets Spec. (y/n)	Measured Time after Cal (m/s)	Passed (Initials)
10	+2%				
Custom Test		Rotor Used:			
	+2%				
	+2%				
	+2%				

#### TEST EQUIPMENT / PRÜFMITTEL

Type	Model	Serial Number	Date of Next Calibration
Stopwatch			
Temperature Meter			
Photo Tachometer			
Safety Tester			

#### SAFETY TEST / SICHERHEITSTEST

Insulations resistance >10 Megaohm	Earth Conductor Resistance <0.2 Ohm	Accessible Current <3500µA	

Revised 14/12/01

Certified By: \_\_\_\_\_ Signature \_\_\_\_\_ Date \_\_\_\_\_ Reviewed By: \_\_\_\_\_ Customer Signature \_\_\_\_\_ Date \_\_\_\_\_ Expiration Date \_\_\_\_\_

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### Labofuge 300

#### Table of Contents – Section 8

TB. No.	Change	Date	Page	Edition
303	Micro Switch of Lid Lock	07/98	7-2/3	01
311	New insulation foil	04/00	7-4	01
320	New main board and lid lock	02/01	7-5	01
329	All versions of the Labofuge 300 have an imbalance switch installed	07/02	7-6	01

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## Technical Bulletin No. 303

### 1. Engineering Change

Introduction of a new micro switch for the lid locking systems of all small centrifuges.

### 2. Reference to Units

Introduction ex works after using up the old stocked micro switches.

Model	Voltage	Cat. No.	Model	Voltage	Cat. No.
Clinifuge	120V	#3538/9	Biofuge 13	230V	#3635/3439
Labo 200	230V	75003630		120V	#3636/7/13/531
	120V	75003631		240V	75003645
	100V	75003597		100V	#3647/3533
	100-120V	75015431		220-240V	75003530
Labo 200S	230V	75003548	Haemofuge	230V	75003520
Labo 200E	220V	75003549		220V	75015358
Medifuge	120V	#3632/11		120V	#3521/22/3612
Medi 200	230V	75003476		100V	75003519
Medi 200S	230V	75003474	Medi 13000	220-240V	75003479
Medi 200E	230V	75003475	Medi 15000	220-240V	75003480
	120V	75003552	Biofuge 15	230V	75003600
Ultrafuge	220V	75003633		220V	75015342
	100-120V	75003636		120V	#3601/4/19
				100V	75003586
Biofuge haemo	230V	75003577	Labo 300	230V	75003230
	120V	75003557		120V	75003231
Biof. fresco	230V	#5510/41	Biofuge pico	230V	#3280/88/97
	120V	#5515/25		120V	#3282/84/89/98
	100V	#5535/6		100V	#3283/96

### 3. Reason for change

The old micro switch with plastic lever is no longer produced. The new micro switch has a metal lever and also improved conditions (closeness) against environmental influences.

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#### 4. Validly

This change is only to perform on units before the above mentioned units, if the micro switch should be defective (e.g. in case of an „OPEN“ or „OP“ message although the lid is closed correctly).

#### 5. Modification

See Service Manual of referring unit: Change of Components for lid lock assembly. Although some Manuals are pointing out that only the complete lid locks are available as one spare part, the new micro switch can be ordered separately, now.

Reinstall the unit accordingly to regulations!

##### Note!

If contrary to expectation a new micro switch should be already defective, so send it back immediately to factory for inspection!

#### 6. Required / Outgoing Parts

Required	Number	Cat. No.	Outgoing/Change	Cat. No.
micro switch per lock (with metal lever)	1	20180557	micro switch (with plastic lever)	20180557

#### 7. Stockage / Obsolete Parts

Old and perfect micro switches can be used up furthering.

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### Technical Bulletin No. 311 / U70

#### 1. Engineering Change

Introduction of a new insulation foil for the Labofuge 300.

#### 2. Reference to Units

This change is effective with the beginning of:

Model	Edition	Cat. No.	Serial No.
Labofuge 300	230V 50/60Hz	75003230	284323
Labofuge 300	120V 50/60Hz	75003231	285898
Labofuge 300	100V 50/60Hz	75003232	285904
Labofuge 300	230V 50/60Hz	75003233	285905

#### 3. Reason for Change

If a lot of liquid comes out into the rotor chamber, it is possible that a part of the liquid comes on the main board. This will prevent by the new insulation foil.

#### 4. Validity

This modification has to carry out if other work (servicing or repair) at this unit is necessary.

#### 5. Modification

Free modification (material) by the Service. Informed customers can order the new insulation foil from Kendro by free of charge.

#### 6. Required / Outgoing Parts

Required	Part No.	Change	Part No.
Insulation foil	20056749	Insulation foil	20056515

#### 7. Stocking / Obsolete Parts

Scrap old insulation foils (PN 20056515) if there are some on the stock.

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### Technical Bulletin No. 320

#### 1. Engineering Change

Introduction of a new main board and latch in Labofuge 300.

#### 2. Reference to Units

This change is effective with the beginning of:

Model	Edition	Cat. No.	Serial No.
Labofuge 300	230V 50/60Hz	75003230	40036367
Labofuge 300	230V 50/60Hz	75003233	40051536

#### 3. Reason for Change

Technical Updating

#### 4. Validity

In case of a failure at PCB or latch system a new PCB or latch **must** be installed.

#### 5. Modification

See replacement description of main board and latch.

#### 6. Required / Outgoing Parts

Required	Part No.	Change	Part No.
PCB 230 V	20150194		
CPU 484	On board		
NV RAM V01	70901073		
Latch 230 V	20901130		

#### 7. Stocking / Obsolete Parts

All units below these serial numbers use the old board and latch further on.

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## Technical Bulletin No. 329 / U82

### 1. Engineering Change

All versions of the Labofuge 300 have an imbalance switch installed

### 2. Reason for Change

Technical updating

### 3. Reference to Units

This change is effective with the beginning of:

Model	Cat. No.	Serial No.
LABOFUGE 300 230V 50/60HZ	75003230	40197991
LABOFUGE® 300 120V 50/60HZ	75003231	40205173
LABOFUGE 300 100V 50/60HZ	75003232	40190446
LABOFUGE 300 230V 50/60HZ SKAND	75003233	40213750
Labofuge 300 100V 50/60Hz SORVALL	75003234	n.n.
Hitachi CT4i 100V 50/60HZ	75003253	n.n.

### 4. Validity

Pay attention to this change on units with higher than above mentioned Serial No's. Note this modification if any work (servicing or repair) at this units is necessary, see new Labofuge 300 Service Manual.

### 5. Modification

This modification cannot perform on units below the above mentioned Serial No's.

### 6. Appendix

New Labofuge 300 Service Manual

# Kendro Service & Support Service Bulletin

<b>Produktlinie</b>	<b>Gerät</b>	<b>Bulletin</b>	<b>Datum</b>
Centrifuges - Heraeus	Labofuge 300	0409-OHA-CEN-349	21 September 2004

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# Kendro Service & Support Service Bulletin

<i>Product Line</i>	<i>Equipment</i>	<i>Bulletin</i>	<i>Issue Date</i>
Centrifuges - Heraeus	Labofuge 300	0409-OHA-CEN-349	21 September 2004

## Change over to new plug system „JST“

**1. Engineering Change:**

Old black plug system (Burndy) becomes obsolete – replacement by new white system (JST)

**2. Reference to Units:**

This change is effective with the beginning of:

Equipment	Version	Cat. No.	Serial
Labofuge 300	230V 50/60Hz	75003230	40235250
Labofuge 300	120V 50/60Hz	75003231	40230962
Labofuge 300	100V 50/60 Hz	75003232	40383659
Labofuge 300	120V 50/60Hz Sorvall	75003233	40237571
Labofuge 300	100V 50/60 Hz	75003234	n.n.

**3. Validity:**

Units **below** the above-mentioned serials require board and latch named **“old”**.

Units **above** given serials require board and latch named **“new”**

**4. Modification:**

Replace main board and latch, see Service Manual Section 5

**5. Required / Outgoing Parts:**

Equipment	Main board “old”	Main board “new”	Latch “old”	Latch “new”
75003230	20150194	20151194 (230V)	20901130	20901143
75003231	20150193	20151193 (120V)	20901131	20901144
75003232	20150165	20151165 (100V)	20901121	20901145
75003233	20150194	20151194 (230V)	20901130	20901143
75003234	20150165	20151165 (100V)	20901121	20901145

**6. Stocking / Obsolete Parts:**

Old PCB 20150193 and 20150194, as well as 20150165 and old latch must be used at older equipment with adapter cables.

<p><b>Author</b></p> <p>Heise</p>	<p><b>Reviewed by</b></p>	
-----------------------------------	---------------------------	---